


PRISM SUPREME HD CUSHION



The Prism Supreme Heavy Duty™ Cushion was designed for users weighing up to 500 lbs and requiring high positioning support. The anatomically designed, firm polyurethane layer, with its aggressive leg troughs and abductor contouring provides superior pelvic stability and leg positioning. A base layer of pressure relieving visco-elastic foam under the IT's (Ischial Tuberosity), allows for excellent immersion for bony prominences.

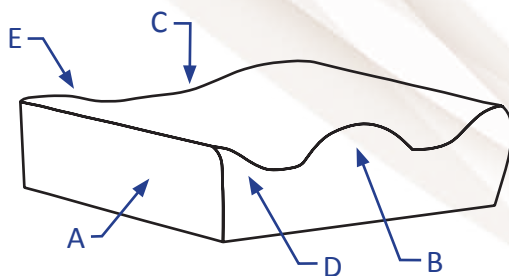

CUSHION THICKNESS

- A- Height at adductor = 6 1/2"
- B- Height at abductor = 5"
- C- Height at rear IT = 4 1/4"
- D- Leg Trough = 4"
- E- Height at rear trochanter = 5 1/4"

Prism Series Cushions		STYLE: Supreme HD	
Stock Sizes – Same Day Shipment			
ADP CODE: SESFM0125			
Part Number	Size W x L	Part Number	Size W x L
SCHD2416	24" x 16"	SCHD2820	28" x 20"
SCHD2418	24" x 18"	SCHD3016	30" x 16"
SCHD2420	24" x 20"	SCHD3018	30" x 18"
SCHD2616	26" x 16"	SCHD3020	30" x 20"
SCHD2618	26" x 18"	SCHD3216	32" x 16"
SCHD2620	26" x 20"	SCHD3218	32" x 18"
SCHD2816	28" x 16"	SCHD3220	32" x 20"
SCHD2818	28" x 18"		


Custom Sizes - Shipment within 7 business days	
ADP CODE: SESFM0125	
Part Number	Size W x L
SCHDCUS	Custom





Product Highlights:

- ✓ High density, contour cut layered foam for a consistent, long lasting product
- ✓ Cane cut in rear to allow for correct placement in wheelchair and proper seat depth
- ✓ Radius cut to match shape of back upholstery
- ✓ Layer of visco foam along the bottom for further pressure relief under IT's (Ischial Tuberosity)
- ✓ Supreme pressure relief
- ✓ 4-way stretch, water resistant inner liner
- ✓ 4-way stretch, breathable outer cover
- ✓ **Product Weight:** 6.8 lbs (26"W x 20"L)
- ✓ **User Weight:** 500 lbs
- ✓ **Warranty:** Two (2) years limited on cushion and 90 days on the cover

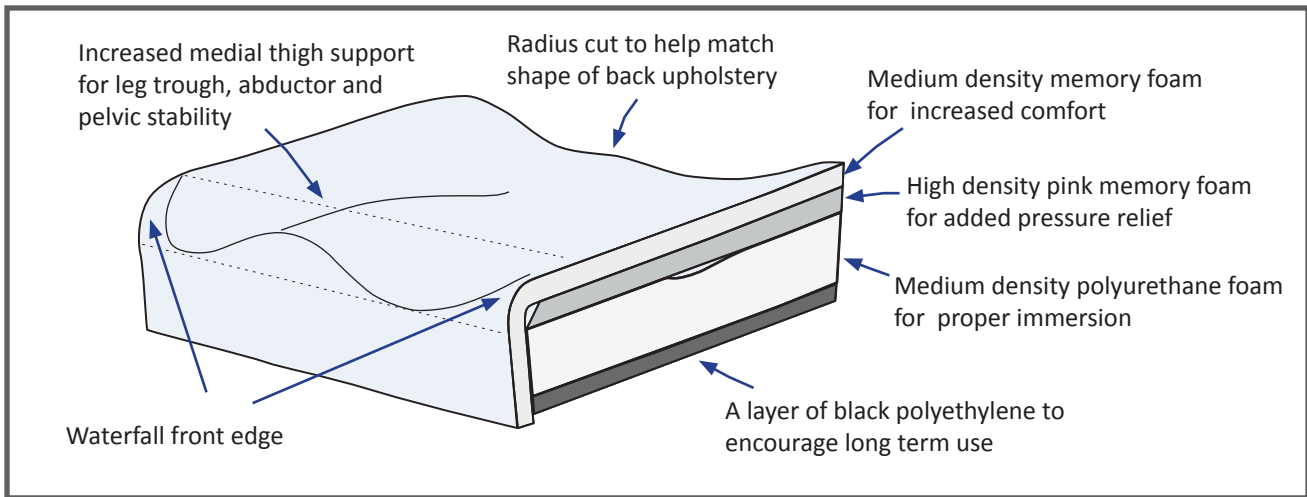

Prism Series Cushions				STYLE: Supreme Heave Duty			
Non-Stock Sizes – Shipment within 3 business days				ADP CODE: SESFM0125			
Part Number	Size W x L	Part Number	Size W x L	Part Number	Size W x L	Part Number	Size W x L
SCHD2422	24" x 22"	SCHD2624	26" x 24"	SCHD2916	29" x 16"	SCHD3118	31" x 18"
SCHD2424	24" x 24"	SCHD2716	27" x 16"	SCHD2918	29" x 18"	SCHD3120	31" x 20"
SCHD2516	25" x 16"	SCHD2718	27" x 18"	SCHD2920	29" x 20"	SCHD3122	31" x 22"
SCHD2518	25" x 18"	SCHD2720	27" x 20"	SCHD2922	29" x 22"	SCHD3124	31" x 24"
SCHD2520	25" x 20"	SCHD2722	27" x 22"	SCHD2924	29" x 24"	SCHD3222	32" x 22"
SCHD2522	25" x 22"	SCHD2724	27" x 24"	SCHD3022	30" x 22"	SCHD3224	32" x 24"
SCHD2524	25" x 24"	SCHD2822	28" x 22"	SCHD3024	30" x 24"		
SCHD2622	26" x 22"	SCHD2824	28" x 24"	SCHD3116	31" x 16"		